

MULTIFUNCTIONAL ICE COOL HAIR REMOVAL INSTRUMENT

Manual



Beauty your life

PRINCIPLE

IPL photon hair removal uses a specific wavelength, the light penetrates the epidermis reaches the straight skin layer, absorbs a lot of photons to generate heat, and is selectively absorbed by the melanin in the hair follicle. The temperature of the hair follicle rises rapidly destroying the hair follicle, and the photons absorbed by the surrounding tissue are not affected. Injury, the growth of hair is delayed or even completely stopped in a short period of time, so as to achieve the purpose of permanent hair removal. Similarly, it can restore the original elasticity of the skin, eliminate or reduce wrinkles, and reduce pores. Improve skin texture, complexion and firm skin.

HAIR REMOVAL

Using the penetration of strong pulses, it directly reaches the hair follicles, causing the hair follicles to generate instant high temperature, causing the hair follicles to be damaged, and the hair growth rate is slowed down to achieve the effect of hair removal.

PHOTO REJUVENATION

For abnormal pigmentation or abnormal skin blood vessels, photorejuvenation phototherapy uses IPL (strong pulse light) to improve the elasticity of the skin and make the skin look younger. This is due to the IPL promoting the synthesis of collagen in the skin and making the skin look fully elastic and shiny.

PDT

In the process of PDT (photodynamic therapy) treatment of acne, IPL is used as a light source to activate substances to shrink the sebaceous glands.

ABNORMAL SKIN AND BLOOD VESSELS

As other laser treatment methods, IPL can also effectively treat vascular structure abnormalities and skin vascular lesions such as flush tubes.

A LOT OF ACNE AND PIMPLES

During the treatment of IPL, a compound porphyrin produced by acne bacteria will appear. When porphyrin absorbs light with a wavelength of about 400 nm, free radicals will be generated, causing damage to acne bacteria. IPL is also effective in treating acne scars.

ABNORMAL SKIN PIGMENTATION

IPL has a significant therapeutic effect on various types of pigmentation abnormalities such as scars, freckles and age spots. Compared with other laser treatments, IPL treatment produces less scabs, so you can apply makeup or wash your face on the day of treatment.

	Laser	IPL
Wavelength	Single wavelength	Many different wavelengths
Treatment	Treatment of 1-2 skin injuries	Treatment of multiple skin injuries
Recovery time	About 1 week	No recovery time
Pain	Local or general anesthesia	No anesthesia
Side effect	Pigmentation	No side effects
Effect	Effective	Without repeated treatment

NO PAIN AT FREEZING POINT

Ice sensation and hair removal go hand in hand, and enjoy the new ice sensation for painless hair removal.

INTELLIGENT RECOGNITION

Adjust the gear according to the state of the skin and hair to remove hair safely and comfortably.

LED INDICATOR

You can enjoy the beauty salon-grade hair removal at home with different lights remind.

PRODUCT DESCRIPTION

Multifunctional hair removal device (hereinafter referred to as hair removal device) is an innovative beauty product. The hair removal instrument has functions of hair removal, skin rejuvenation and acne removal, and adopts intelligent pulse light irradiation technology, which integrates multiple applications of optical technology in the field of beauty and health. It is one of the most professional and practical products in this field.

PLEASE READ THE INSTRUCTIONS CAREFULLY BEFORE USE

The design of OSENYUAN hair removal instrument is always based on the customer's comfort and safe use. It is manufactured under the highest quality standards. Before using the product, please read and understand all instructions and precautions.

APPLICABLE PEOPLE

Hair removal device is an over-the-counter device suitable for men and women to remove excess hair.

APPLICABLE SKIN COLOR

Black skin is prohibited.

Please use the brown skin from the smallest level, and adjust it to the corresponding field according to your skin's ability to adapt to the pain.

SAFETY PRECAUTIONS

CONTRAINDICATIONS

The use is prohibited in the following situations:

Children and elderly people, local or systemic inflammation; people with photosensitive skin and people who are using photosensitive drugs, users of afferent medical electronic devices such as heart pacemakers, artificial heart and lung and other life-sustaining devices may cause electronic medical treatment Incorrect operation of the device.

MATTERS NEEDING ATTENTION

1. Pregnant women, under 11 years old, black moles, tattoos, etc. are not recommended. It is not recommended to use sunbeds or sunbeds and sunless tanning products. If you are not sure whether your skin is tanned or sunburned, we recommend waiting a day before using it.
2. Do not use tattoos, dark brown and black spots, these parts may absorb a lot of light and cause skin discomfort.
3. It is forbidden to use around the ears, neck, nipples, genitals or anus. These parts of the skin are sensitive and vulnerable to damage.
4. Do not stay in the same position for too long, or emit light multiple times in one position, which may generate too much heat and make the skin uncomfortable. If the skin is red and swollen after use and lasts for more than 24 hours, please suspend the use, please contact the doctor or consult customer service.


SAFETY WARNING







1. There is high voltage in the hair removal instrument, and it is strictly forbidden to disassemble it without permission.
2. When using the epilator, please do not cover the vents.
3. The hair removal device must not get wet or water.
4. It is forbidden to use the epilator when the plug or power cord is damaged.
5. It is forbidden to operate the epilator in flammable and explosive environments.

SAFE OPERATION SPECIFICATIONS

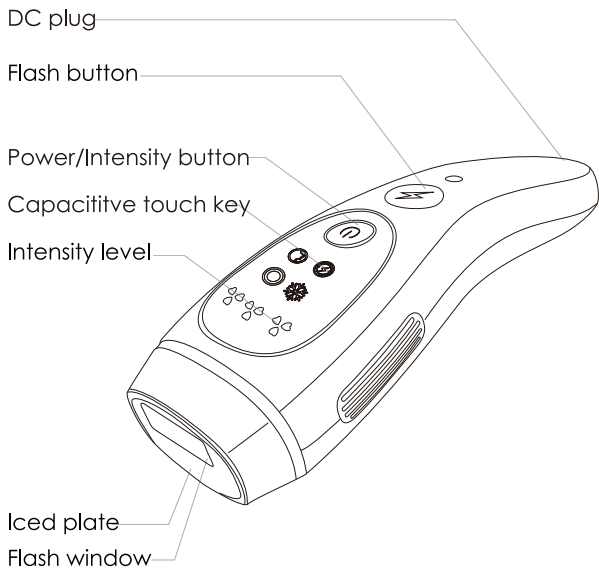1. During flash output, it is forbidden to point the bald head at the fragile parts such as the eyes. Wear goggles when using.
2. Before using this product, please remove all metal objects carried on the body, such as ring bracelets, etc., so as not to cause equipment failure.
3. Please operate in a place away from children. After use, it should be placed in a place where children can't reach it.

BUTTON FUNCTION INTRODUCTION

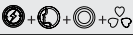

LED INDICATOR	FEATURES	PARAMETER
	Hair removal	999999 flashes
	Skin rejuvenation	
	Acne removal	
	ICE Cool	
	Automatic mode	Flash
	Intensity display	1-9
/	Remaining times	/

SKIN SENSITIVITY TEST

If you are using this product for the first time, perform skin sensitivity test according to the following requirements:

- ① Turn on the computer according to the following steps.
- ② Select the gear 1st gear or 2nd gear.
- ③ Operate the selected skin area according to the operation instructions and precautions.
- ④ After 30 minutes of skin test, check the skin reaction.
- ⑤ If the skin of the test site is normal, you can continue to use it. If there is abnormality, please stop using it.

STEPS FOR USAGE

Power on: Plug it in, all lights blink once (about a second)

After all buttons flashing is completed, the power button still flashes (representing the state that can be turned on)

1. Start: Press the power button once,

2. After power on, the hair removal mode + ice button + power button is always light on + flash button is always light on

3. After touching the skin, the level light flashes, the flash button goes out after each flash, and then the internal energy storage capacitor is fully charged and then lights on again

4. Long press the Flash button

5. Hair removal mode + level light flashes (PS), the power button is always light on; the flash button is off after each flash, and then the internal energy storage capacitor is fully charged and then always light on again, totally nine levels.

6. Acne removal mode: replace the filter head (560nm filter), touch the acne removal mode light, and the acne removal mode light flashes at the same time, press the power button to shift, and one button corresponds to an increase of a gear, a total of nine Gear,

7. Facial beauty mode: replace the filter head (600nm filters), touch the face beauty mode light, and the face beauty mode light flashes at the same time, press the power button to shift, and one button corresponds to an increase of a gear, a total of nine Gear.

The gear display light flashes after the touch is sensed to the skin; press the flash button once to flash once; (long press cannot flash continuously when the skin is sensed)

Automatic mode - (long press the flash button to enter automatic mode - need to operate without touching the skin, and the hair removal mode light is always on)

8. Hair removal function: press the power button to shift, and one button corresponds to increase a gear, a total of nine gears, 510nm filter;

9. Acne removal mode: replace the filter head (560nm filter), touch the acne removal mode light, and at the same time the acne removal mode light is always on, press the power button to shift, and one button corresponds to increase a gear, a total of nine gears.

10. Facial beauty mode: replace the filter head (600nm filters), touch the face beauty mode, and the face beauty mode light is always on, press the power button to shift, and one button corresponds to an increase of a gear, a total of nine Gear.

Ice compress function: Press the function ice compress button to start the ice compress function; the ice compress function in acne removal mode and facial beauty mode is disabled;

PS: The level button flashes when it touches the skin, and is always light on when it does not touch the skin

11. Power off: long press the power button to power off, except for the power button is flashing, all other lights are off

12. No sound to turn off immediately: Unplug the plug, the power button is light off

NOTE

1. After removing hair, apply aloe vera gel or similar moisturizing cream to the hair removal area. Do not use harsh cosmetics immediately.
2. Be sure to clean the treatment area before hair removal. It is recommended to wash the skin with cold water. Do not touch the water within 2 hours after the operation. Do not sunbathe, warm water bath, steam bath, etc, within 48 hours.
3. Avoid exposing the hair removal site to strong sunlight within 2 hours after hair removal and apply sunscreen before outdoor activities. It is recommended to operate at night.
4. Do not swim in the disinfected swimming pool within 48 hours after hair removal.
5. Do not press the trigger button continuously, just press it once. It is recommended to irradiate each part 2-3 times. (The area of the skin where the entire lamp head is in close contact is one place).
6. The continuous working time of the machine should not exceed 25 minutes. If the hair removal area is relatively large, it can be divided into several times (the power should be turned off after 10 minutes, and you can continue to use it after 10 minutes of rest)

USAGE TIME

The hair has a certain growth cycle, and short-term use cannot complete permanent hair removal. The entire treatment process needs to be carried out for 2-6 months. Only long-term adherence can see the effect.

First month (use once every 3-5 days)

After the irradiation, the hair follicles began to shrink, and the hair gradually became soft and thin.

Second month (use once a week)

Continuously irradiated hair has slowed growth, and hair invention has been significantly reduced.

Third month (use once every 15 days)

The hair roots of the hair follicles are irradiated according to the cycle, and no dark spots are left, and the skin is silky.

Fourth month (use once every 30 days)

96% of hairs stop growing, achieve permanent hair removal, and have perfect skin.

Note: The hair removal instrument is used to remove hair for a long time.

Some parts do not require subsequent maintenance, and some parts require subsequent maintenance at different times of 6 months, 9 months, and 12 months.

Function switching

- ① Follow the steps of using the epilator to the LED indicator interface
- ② Ice compress function: short press "switch on/off button", the snowflake icon appears on the display, and then enter the ice compress function
- ③ Hair removal & skin rejuvenation & acne removal function: single press "function switch button", the display will switch at "HR" "SC" RA

Suggestions for the three functions of hair removal/skin rejuvenation/acne removal

Each part can only be irradiated with one function photon in a day. If the same part needs to be irradiated with other functions, it needs to be operated every other day.

Precise hair removal

The machine is additionally equipped with round and square small lamp heads, which can accurately remove hair from small parts such as fingers. When using, just place the small lamp head on the flash outlet of the host, it will automatically adsorb on the host and can be used directly Off

SPECIFICATIONS

Light source	IPL	power	100-240 V, 50-60Hz, 48W
Light area	3cm ²	Panel	Button + touch button
Host size	212*63*67mm	Function	HR/SC/RA
Total weight	Standard as package list	work	Temperature -20-60° relative humidity 20-90%
Machine weight	270g	Warranty	One year warranty

PACKING LIST

Host	1
Shaving knife	1
Instruction	1
Glasses	1
Packing box	1

ADVERSE REACTIONS AND TREATMENT


All beauty methods may carry certain risks, and this product may cause discomfort when used improperly. The problems and suggestions are as follows

Tingling sensation	Everyone has different pains because of their different constitutions. If there is a tingling sensation, please lower the gear
Slight burning sensation	The normal phenomenon can disappear by itself. You can use aloe vera gel or a towel to cool down.
Burnt scent	Normal phenomenon, be sure to shave hair before using the epilator

SIMPLE TROUBLESHOOTING

Phenomenon	Possible reason	Processing method
Do not show light	1. The switch button is not pressed into place 2. Poor power contact	Connect the power cord
Not boot No flash	1. No contact with the lamp head and vertical skin 2. The number of lamp head flashes has been used up	1. Close to the skin vertically 2. Press when the ready indicator flashes
Automatic shut-down	The system is set to not flash for 5 minutes and automatically shut down	Switch power

BODY PARTS

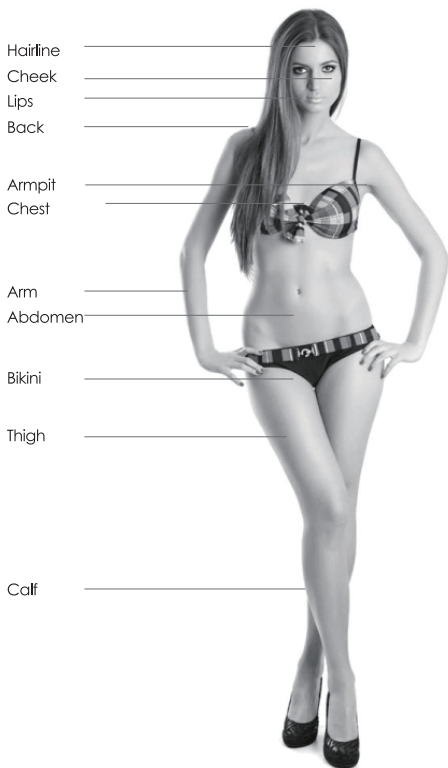

Hairline

Cheek

Lips

Back

Armpit

Chest

Arm

Abdomen

Bikini

Thigh

Calf

FREQUENTLY ASKED QUESTIONS

Q: Should I shave before each operation?

A: No, just scrape once when using it for the first time. After that, you can wait for the hair to grow to about 5 mm. If the effect of use is affected, just scrape it again.

Q: Why didn't my hair take off immediately after irradiation, or did it change after irradiation?

A: This is a permanent hair removal device. You need to use a course of treatment before you can permanently remove hair. At the beginning, it can only inhibit hair growth and slow down the hair growth rate. Under continuous action, the hair can be put into a dormant state to achieve permanent hair removal effect.

Q: Why did I use it so few times that it didn't work as well as others, or no effect?

A: Everyone's physique is different, the thickness of the hair is different, the parts used are different, and the effect cycle is naturally different. The use cycle is an average value, and the specific time may be slightly different.

Q: What gear is appropriate?

A: There is no unified standard, depending on everyone's skin tolerance (the darker the skin color, the lower the gear used). Slowly increase the use from the lower gear, the higher the gear, the greater the energy, the more obvious the effect of use. The treatment time will be shortened accordingly. On the contrary, if the low gear is used, the treatment time will be correspondingly extended. It is recommended that each part of the skin be tested for 1-3 flashes of light on the skin first. After the skin is free of discomfort on the second day, use it in a large area.

Q: Will it hurt?

A: During the hair removal process using the hair removal instrument, light waves will encounter heat melanin and generate heat energy, so there will be some subtle pain. Everyone's skin sensitivity is different, and the ability to withstand pain is different, so every Personal feelings are different, it is recommended to increase the gear slowly from the low gear according to your actual situation.

Q: How to shave? Which hair should be shaved?

Answer: Hair like armpits, bikinis, legs and other parts should be shaved first and then used. (Recommendation: Long shave the hair you need to remove in advance to less than 1mm.) Before shaving, it is recommended to use the following two This method makes shaving easier

a: Lubricate the shaved area with shaving cream or mild shower gel

b: Use a hot towel to apply to the shaved area for a few minutes, and gently shave the hair immediately after application; this will make the operation easier.

Q: Why do I have a scorched smell after use?

A: The hair is not completely shaved during use. Melanin will interact with light waves and produce heat energy, so it will have the smell of burning hair. This is normal.

Q: I want to remove my hair quickly. Can I use it every day?

A: It is not recommended to use it every day. There is a cycle of hair growth. You can use it according to the normal use cycle. You can use it once every 3-5 days during the first month. After the hair growth rate slows down, the number of times will gradually decrease.

Thank you very much for purchasing
IPL hair removal
Hope our products and services will
make you satisfied

